


# NEUATION

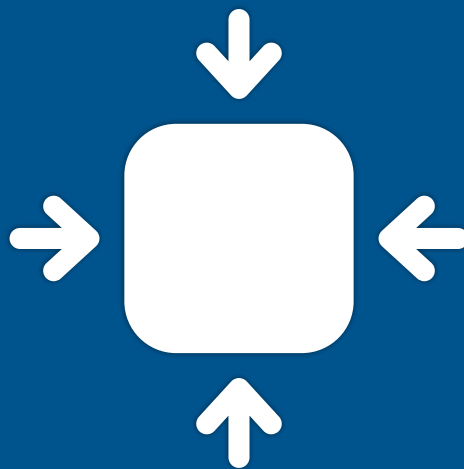

COMPACT DEVICES  
THAT OFFER  
IMMEDIATE VORTEXING

iSwix Jr. NXT and iSwix Jr. VT NXT

Mini Vortex Mixer

## iSwix Jr. NXT & iSwix Jr. VT NXT

### Mini Vortex Mixer

Mini vortex mixer iSwix Jr. NXT and iSwix Jr. VT NXT are compact devices that offer immediate vortexing for a variety of tube sizes. Comes with a simple touch start operation enabling exceptional mixing performance with one hand.

Available in two variants: Fixed speed and variable speed with a maximum **RPM of 4200**. Its powerful motor provides a smooth, quiet and efficient mixing without any walking issue



## Product Features



Touch Operation-can handle plastic tubes to up to 30mm in diameter




Extremely small footprint



4.5mm orbital motion

Suction cups for proper grip

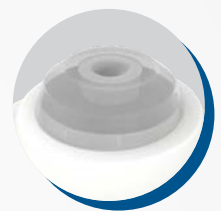

Low voltage power supply



Suspension design for longer switch life



Variable speed for gentle to vigorous mixing (Variable model)



# Specifications

iSwix Jr.	Fixed	Variable
Motor Type	DC	
Type Of Movement	Orbital	
Maximum Speed	4200RPM	2000-4200RPM
Shaker Diameter	4.5 mm	
Tube Diameter	Up to 30 mm	
Speed Accuracy	± 250 RPM	
Operation Type	Touch	
Protection Class	IP 43	
Permissible Ambient Temperature	4°C to 45 °C	
Weight	0.75 Kg	
Permissible Relative Moisture	<80%	
Noise Level	<50 dB	
Dimensions (L X W X H)	120 x 120 x 103mm	
Input Voltage	100-240VAC 50/60Hz	
Power Consumption	52W	

# Ordering Information

iSwix Jr. NXT/iSwix Jr.VT NXT		
Product Name	Cat#	Detail
iSwix Jr NXT.	Jv4210 EU/US/UK/AU	Fix Speed Mini vortex mixer plastic body - Max RPM 4200
Iswix Jr VT NXT	Jv5210 EU/US/UK/AU	Variable speed Mini vortex mixer plastic body - Max RPM 4200

# NEUATION

www.neuation.com

AN ISO 9001 : 2015 & ISO 13485 : 2016 CERTIFIED MFG. COMPANY

\*We also execute OEM projects. Specifications can be changed without notice for quality improvement.