




Cirrus™ 2

A range of nebulisers and kits

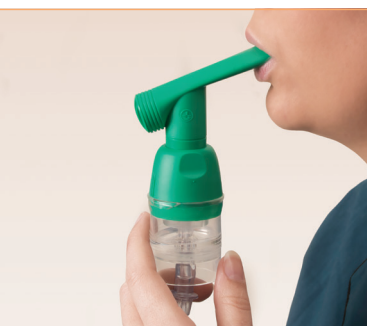

Oxygen and Aerosol Therapy ▪ Nebuliser Therapy



Cirrus™2 nebuliser and nebuliser mask kits

Based on the well-established Cirrus nebuliser the Cirrus2 incorporates a number of improved features. These include a lower residual volume to reduce drug wastage and a graduated drug cup, which provides guidance on how much drug is in the nebuliser.

The cup now has a quick-release fitting which enables the nebuliser to be attached, as well as detached, with a quarter turn.


This speeds up the refilling process and provides a more

▶ Video available at www.intersurgical.com

reliable seal between the cup and the nebuliser body.

The Cirrus2 nebuliser is designed to deliver drugs for tracheobronchial deposition. At a driving gas flow of 8 L/min the Mass Median Diameter (MMD) is 3.3 microns.

The Cirrus2 Adult mask kit incorporates the Intersurgical EcoLite™ aerosol mask, this non-PVC mask reduces the environmental impact of this product compared with a conventional mask kit by 28%.

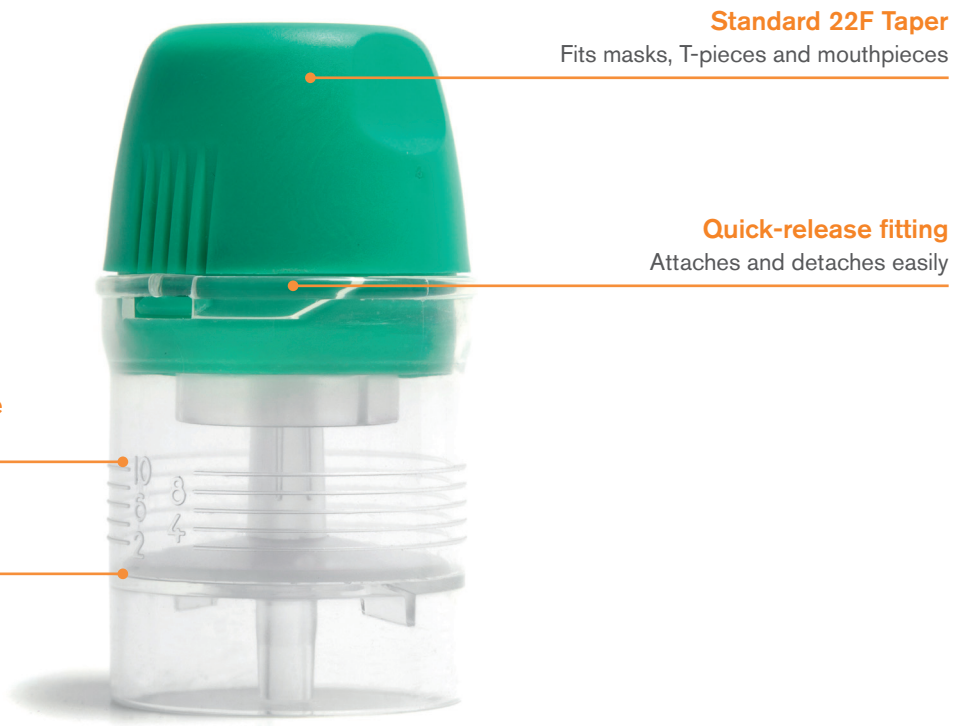

Standard 22F Taper

Fits masks, T-pieces and mouthpieces

Quick-release fitting

Attaches and detaches easily

Graduated nebuliser cup guide

Faster preparation

Low residual volume

Reduced drug waste

Self-sealing T-piece

The Cirrus2 is now available with a self-sealing T-piece, which enables the nebuliser to be attached (or detached) from the ventilator without interfering with patient ventilation or inadvertently activating any alarms.

▶ Video available at www.intersurgical.com





Mask kit with the new Intersurgical EcoLite™ aerosol mask

The Intersurgical EcoLite aerosol mask has been designed as part of our ongoing focus on sustainable development, which is an integral part of our ethos.

The Eco range, unlike traditional masks, contains no PVC which provides an option with a significantly lower impact on the environment.

Improved patient comfort – Lightweight and soft, patient comfort is also paramount in the design.

By utilising the latest technology we have combined two materials: the first forming the clear rigid shell of the mask, this provides a strong lightweight body with good visibility, whilst the second softer material forms the seal, providing a comfortable fit for a wide variety of face shapes.




Incurved nose seal

Conforms to different nose shapes designed to prevent aerosolised drug patient's eyes

No metal nose clip

MRI compatible

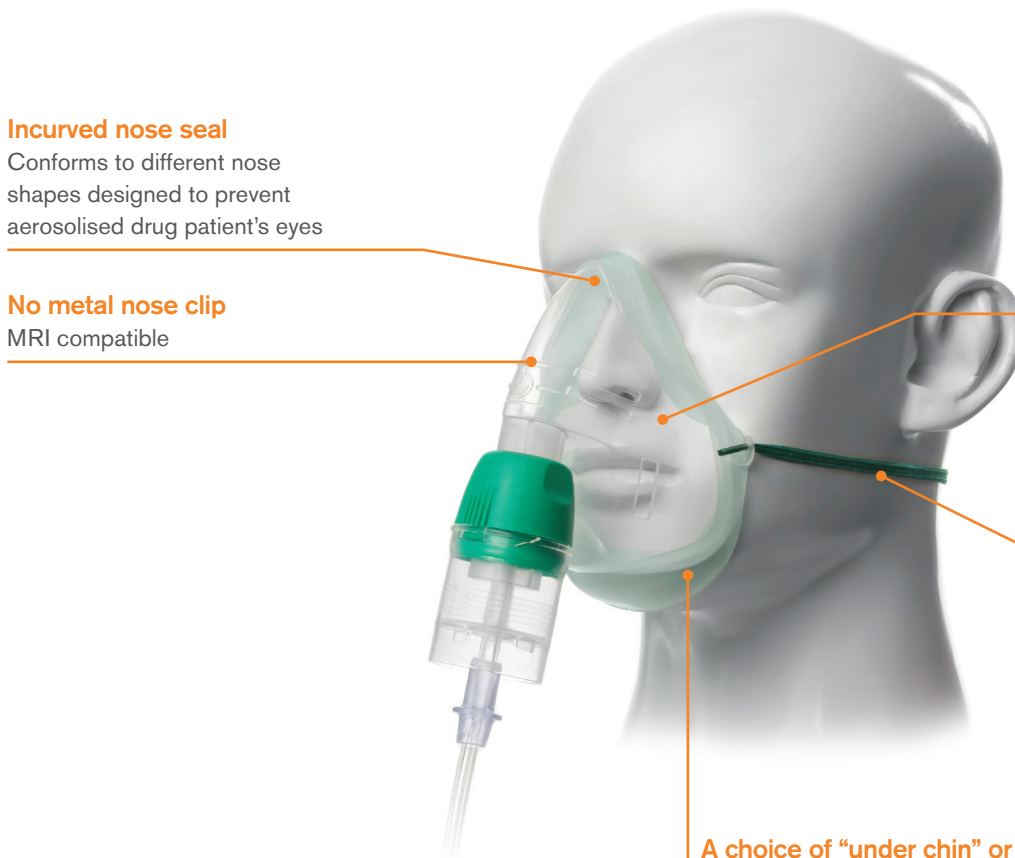

Soft face seals

Increased patient comfort

Elastic can be positioned under or over the ears

Below ear position eliminates trauma to top of ears

A choice of "under chin" or "on chin" positions

Provides a better fit on a wider range of patient face shapes

▶ Video available at www.intersurgical.com





Comfortable for the environment

Various groups have questioned the use of PVC in medical products and it's impact on the environment. Intersurgical's aim is to reduce the environmental impact of its products and processes. The utilisation of new materials and the latest manufacturing technology has resulted in the elimination of PVC from the mask shell resulting in a reduced environmental impact.



www.intersurgical.com/info/eco

Cirrus™2 nebuliser and nebuliser mask kits

Code	Description	Tube length	Box Qty.
1402000	Cirrus2 nebuliser		75
1453015	Cirrus2 nebuliser, adult, Intersurgical EcoLite™ mask kit and tube	2.1m	30
1453010	Cirrus2 nebuliser, adult, Eco mask kit and non-PVC tube	2.1m	30
1454015	Cirrus2 nebuliser, paediatric, Intersurgical EcoLite mask kit and tube	2.1m	36
1455000	Cirrus2 nebuliser, one-piece mouthpiece kit and tube	2.1m	35
1464000	Cirrus2 nebuliser, universal mouthpiece T-kit		40
1464004	Cirrus2 nebuliser, universal mouthpiece T-kit and filter		20
1464021	Cirrus2 nebuliser, mouthpiece kit, filter and tube	2.1m	20
2507002	Cirrus2 neonatal nebuliser bagging kit		18
2605000	Cirrus2 nebuliser breathing system T-kit, 22mm		40
2606000	Cirrus2 nebuliser breathing system T-kit, 15mm		35
2608000	Cirrus2 nebuliser breathing system T-kit, 10mm		45
1958000	Cirrus2 nebuliser anti-pollution control kit		20

Self-sealing T-piece

Code	Description	Box Qty.
2609000	Cirrus2 nebuliser self-sealing T-piece kit, 22mm	35
1814000	Self-sealing in-line T-piece 22M-22F, 22M	40

Make an enquiry

Lower environmental impact product



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 9656 300 F: +44 (0)118 9656 356 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 9001:2008, ISO 13485:2003 and ISO 14001:2004

Please think before you print
Save energy and paper.
If you must print this information sheet please print it double sided.

IS10.23 • Issue 10 09.17